使用說明書
INSTRUCTIONS


HITACHI
Inspire the Next

日立冰肌離子清潔及保濕器 (交流式)
HITACHI Ion Cleansing with Moisturizing and Cooling Support Device

CM-N10000BF

只作家居用途
For household use only

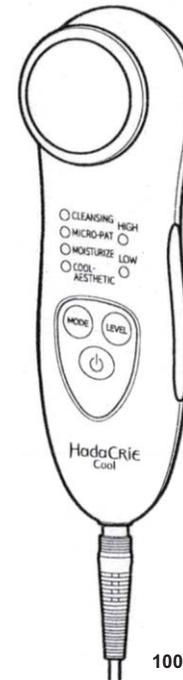100V-240V AC

日本原装
Made in Japan

Hitachi (Hong Kong) Ltd.



18th Floor, Ever Gain Centre, 28 On Muk Street,
Shatin, N.T., Hong Kong

Hitachi Maxell, Ltd.

4680 Ikata, Fukuchi-Machi, Tagawa-Gun, Fukuoka,
822-1296 Japan

English Page 1~12

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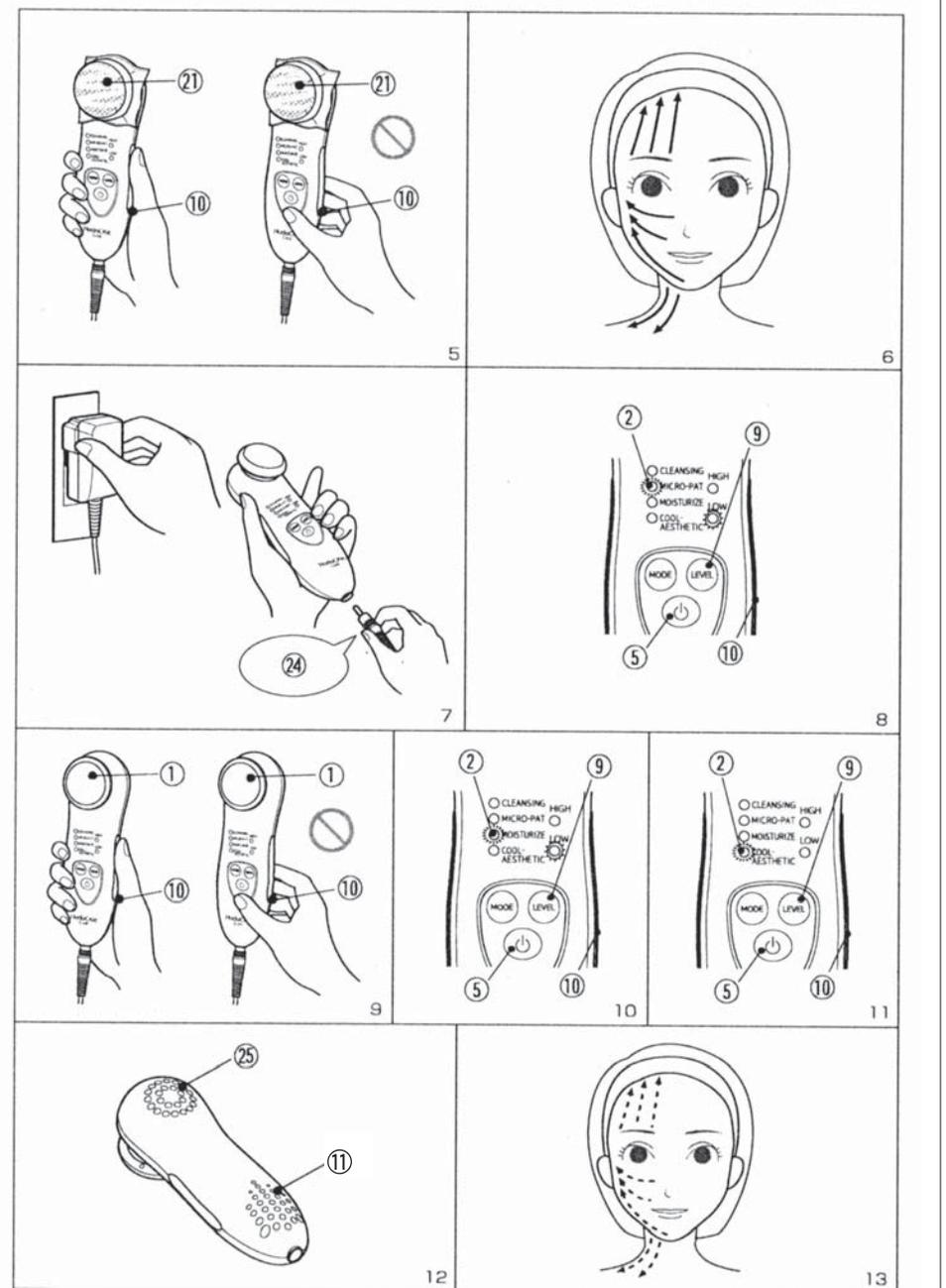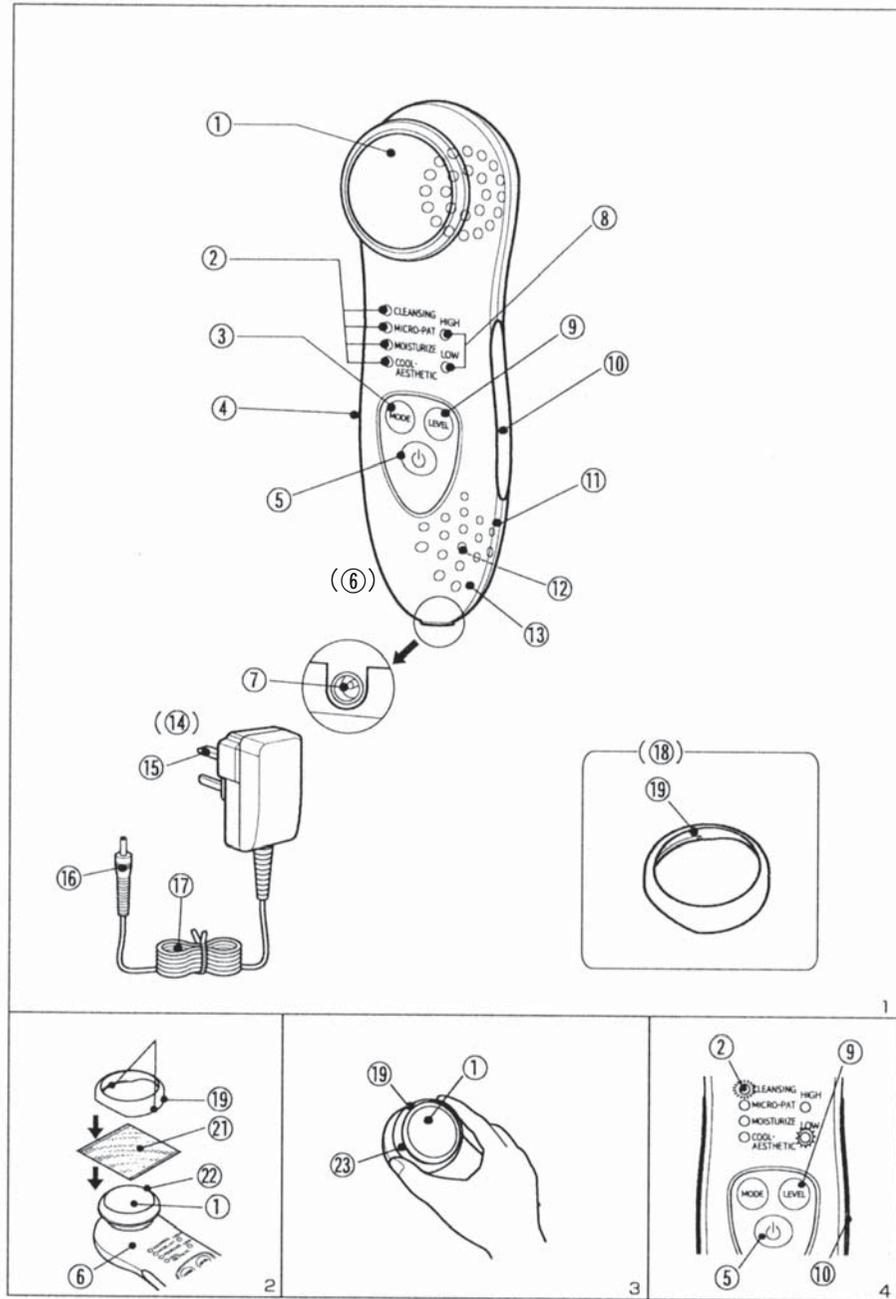English

Safety precautions

- Precautions described here are intended for safe and correct use of this product and are to prevent injury and damage to you and others before they occur.
- Precautions are classified as follows.

Make sure to observe the precautions because they are all important in relation to safety.

 DANGER	Indicates an imminently dangerous situation in which improper handling may result in death or serious injury to a user.
 WARNING	Indicates a dangerous situation in which improper handling could result in death or serious injury to a user.
 CAUTION	Indicates a dangerous situation in which improper handling could result in slight injury to a user or damage to property (house, household effects, etc.).

<Examples of pictorial indication>

	The  symbol indicates "Prohibited" (Do not do this).  "Disassembly prohibited"  "Use prohibited in wet areas"  "Use with wet hands prohibited"
	The  symbol indicates "Compulsory" (Make sure to do this).  "General action"  "Unplug"

DANGER

	Persons who use the following medical electronic devices should not use this product. (This product may cause the medical electronic devices to malfunction.) - Self-contained medical electronic devices such as a pacemaker - Wearable medical electronic devices such as an electrocardiograph
---	---

WARNING

	Do not use this product by other methods which are not described in this instruction manual. (Cause of accident, trouble, or failure)
	Use this product only for face and neck. (Cause of an accident or trouble)
	The following persons should be avoid to use this product. (Cause of an accident or trouble) - Persons with heart disease - Persons with an acute (pain) illness - Persons with fever - Persons with a blood pressure disorder - Persons with an infectious disease - Persons with a tubercular disease - Persons with a malignant tumor - Persons with facial neuralgia - Persons who are pregnant - Persons who are menstruating - Persons with atopic dermatitis - Persons with an allergic predisposition - Persons with particularly sensitive skin - Persons with hemophilia - Persons under dental care - Persons whose blood vessels have expanded due to long time use of steroid hormones or liver function disorder - Persons who are under the influence of alcohol or medication, or who are very tired

WARNING

	Do not use this product on the following regions. (Cause of accident or trouble) - Regions where orthopedic treatment was given - Eyeballs and upper eyelids - Wound sites - Facial melanoderma - Mucosal regions, such as inside the mouth - Regions irritated due to pimples or skin rash - Regions which suffer from sensation of pain or perception disorder - Regions where metal, plastic or silicone etc., is embedded under skin - Regions where there is itching, sensation of warmth, or liver spots due to physical stimulus etc. - Regions where skin irritation is caused by cosmetics
	This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
	When feeling pain or disorder during use of this product, or sensing anything wrong with your skin after use, stop using this product and consult a dermatologist. (Cause of skin trouble etc.)
	If any abnormality occurs during use, immediately pull out the power plug of the power adaptor from the outlet. (Cause of accident or fire)
	Use only 100-240V AC for power supply. (Cause of fire or electric shock)
	Do not use if the power cord or the power plug is damaged or its insertion to the outlet is loose. (If the power cord breaks, a fire or electric shock may occur.) The power cord cannot be unnecessarily bended, pulled, twisted, bound, placed a heavy object on it, pinched, processed, or wound around the power adaptor. (If the power cord damages fire or electric shock may occur.)
	If any unusual sound or abnormality occurs in the main body, immediately stop using this product and contact the distributor for repair together with the power adaptor. (Cause of accident or trouble)
	Persons who use this product for the first time should use it while observing the condition of their skin. (Cause of skin trouble or poor physical condition) Use the COOL-AESTHETIC mode at an ambient temperature in the range of 10 to 40°C. (Cause of skin damage as the head is too cold if the device is used at a low temperature)
	Do not insert or remove the power plug in/from the outlet with wet hands. (Cause of failure or short circuit)
	Do not wash or splash water on this product. (Cause of failure or short circuit)
	Do not use or keep this product in places easily splashed with water or in damp places, such as in a bathroom. (Cause of failure or short circuit)

⚠ WARNING

Do not alter disassemble or repair by unauthorized technician.
(Cause of fire, electric shock, injury, or failure)

⚠ CAUTION

Check the head for any damage or deformation before use.
(Cause of skin damage)

Do not use this product together with other beauty treatment devices at the same time. (Cause of failure or trouble)

Do not use this product if there is some injury or wound to your skin.
(Cause of skin trouble or pain)

Do not drop this product or give it a strong impact.
(Cause of injury or failure)

Do not wipe with benzene, thinner or alcohol, and do not spray insecticide.
(Cause of cracking, electric shock, or ignition)

Use only the attached power adaptor. And, do not use any other product with the attached power adaptor. (Possibility of electric shock or ignition)

Remove accessories such as earrings and necklaces, and contact lenses before use. (Cause of accident or trouble)

If there is any severe skin disorder or irregularity on your skin surface, or due to downy hair, you may feel strong irritation in rare instances. In this case, avoid using this product on that region. (Cause of accident or trouble)

Do not keep the cotton puff or the head on one part of your skin but always keep it in motion. (Cause of skin trouble or poor physical condition)

Do not hold the cord to pull out the power plug and device plug.
(Cause of disconnection, electric shock, short circuit, or ignition)

Make sure to pull out the power plug from the outlet while not in use.
(Cause of electric shock or fire by deterioration)

Names of parts and their functions (Fig. 1)

- ① Head (titanium material)
Set a cotton puff and apply it to the skin.
This section vibrates during use in modes except "COOL-AESTHETIC" mode.
- ② Mode indicator lamp (orange)
- ③ "MODE" selector button
Every time this button is pressed, the mode is switched.
- ④ Side plate

Mode	During preparation	During use
CLEANSING	Orange flashes	Orange lights
MICRO-PAT		
MOISTURIZE		
COOL-AESTHETIC	Blue flashes	Blue lights

- ⑤ Power switch ("⏻")
Press this switch to turn on the "LOW" level indicator lamp lights.
Press this switch again to turn off.
Equipped with an auto-off function, it will automatically turn "OFF" if the device is not used for about 30 seconds in an "ON" state.
- ⑥ Body
- ⑦ Body socket
- ⑧ Level indicator lamps
"HIGH" (red) : "LOW" (green)
When the power switch is turned "ON," "LOW" lights. It does not light in the "COOL-AESTHETIC" mode.
- ⑨ High/low selector button
Every time this button is pressed, the power level is switched between "HIGH" and "LOW." The vibration power does not change.
- ⑩ Side plate
Touch this part and apply the head to your skin, and the device will start to operate in modes except "COOL-AESTHETIC" mode.
- ⑪ Exhaust hole (back)
- ⑫ Cooling fan (built-in)
This operates in the "COOL-AESTHETIC" mode.
- ⑬ Buzzer (selecting sound/ready sound/warning sound/ending sound)
Every time the switch/button is pressed, a selecting sound can be heard.
At ready in the "COOL-AESTHETIC" mode, a ready sound can be heard. If a warning sound can be heard during use, it is a signal to indicate a device failure.
When each mode ends, an ending sound can be heard.
- ⑭ Power adaptor Model: KH-186B
- ⑮ Power plug
Insert in an outlet to the roots.
- ⑯ Device plug
- ⑰ Cord
- ⑱ Accessories
- ⑲ Rings (two)
Rings for fixing a cotton puff.

Before use

⚠ WARNING	When sensing anything unusual with your skin or body during use, immediately stop using this product and consult a dermatologist. (Cause of an accident or poor physical condition)
⚠ CAUTION	Remove accessories and contact lenses. (Cause of skin trouble or poor physical condition)

1. About mode selection

Step	Mode	Cotton puff and cosmetics
1	"CLEANSING" (To remove fine dirt remaining in pores)	Use a new cotton puff saturated with your regular toner.
2	"MICRO-PAT" (To allow moisturizing ingredients to easily permeate your skin)	Apply your regular toner to your skin for use. No cotton puff is used.
3	"MOISTURIZE" (To give moisturizing care by making moisturizing ingredients permeate your skin)	Use a new cotton puff saturated with your regular toner.
4	"COOL-AESTHETIC" (To tighten the pores of your skin by the cooling effect)	Apply the head section directly to your skin for use. No cotton puff or toner used.

How to use

- Make sure to use the attached power adaptor on an AC power supply of 100 – 240V. (This product is not rechargeable.)

⚠ WARNING	Do not use this product around eyes and ears. (Cause of an accident or poor physical condition)
	If you sense anything wrong with your skin or body during use, immediately stop using this product and consult a dermatologist. (Cause of an accident or poor health condition)

Step 1 "CLEANSING"

This mode is to remove fine dirt remaining in pores which cannot be removed by simple washing and make your skin suitable for "MICRO-PAT".

1. After washing your face, saturate the cotton puff attached to the head with toner. Use your regular toner or wiping lotion with an exfoliating effect.

(Do not use a milky lotion or gel.)

- Saturate the cotton puff with the toner until it seeps through the cotton puff to the rear side.

Some toners (made only of mineral spring water or water, etc.) may not show effects sufficiently.

2. Connect the power adaptor to the main body and outlet, and press the power switch. (Fig. 4)
 - The "LOW" level indicator lamp will light, and the "CLEANSING" mode indicator lamp will flash. If another mode indicator lamp flashes, press the mode selector button to select "CLEANSING."
3. Select the power level ("HIGH" or "LOW") with the level selector button. Select the power level ("HIGH" or "LOW") with the level selector button according to your skin condition.
 - Every time the level selector button is pressed, the potential power level is switched between "HIGH" and "LOW". The vibration power does not change.
 - The previous power level ("HIGH" or "LOW") is stored by the memory function.

Use at the "LOW" level is recommended for people who use this product for the first time and who feel something unusual.

4. Hold the main body by hand touching the side plate. (Fig. 5)
 - Make sure to hold the body with your hand touching either of the right and left side plates.
5. Apply the cotton puff to your skin and move the product slowly. (Fig. 6)
 - The mode indicator lamp will light, and the product will start to vibrate.
 - Use this product within approximate 5 minutes.

After 5 minutes have elapsed, the product will stop automatically.

CAUTION

The device will not vibrate unless you apply the cotton puff to your skin while touching the side plate. (Skin sensing function)
If the head is not vibrating, the product is not working.
If you leave the device as it is for about 30 seconds without vibration, the power will turn off automatically. (Auto-off function)

- Lightly apply and slowly move as if sliding from the center to the outside of your face.

⚠ CAUTION Do not keep the cotton puff on the same part of your skin, but move it around. (Cause of skin trouble or poor health condition)

- If the moisture in the cotton puff becomes insufficient or if it is hard to slide, add more toner.
 - Do not use this device on the same part of your skin for more than 10 seconds.
 - If you feel slight current stimulation (a tingling sensation) during use, press the level selector button to set the "LOW" level.
If you still feel a "tingling sensation," add more toner.
6. An automatic stop is made when around around 5 minutes have elapsed. (Auto-off function)
 - You will know the automatic stop with the ending sound.
 7. After use, make sure to pull out the power adaptor from the outlet. (Fig. 7)
 - ④ Pull out holding the plug portion.
For continued use of the next mode, keep the power adaptor connected.
 8. Remove the ring and the cotton puff.
 - Wipe off the surface of the head with tissue, etc.
 - Do not use the cotton puff again once it is used because its reuse will cause skin trouble.

Step 2 "MICRO-PAT"

This mode is to effectively stimulate and loosen up your skin to allow moisturizing ingredients to easily permeate your skin.

1. Evenly spread your regular toner over your face.
 - Wipe the toner off your hands.
2. Check that the power adaptor is connected to the main body, and press the power switch. (Fig. 8)
 - The "LOW" level indicator lamp will light, and the "MICRO-PAT" mode indicator lamp will flash. If another mode indicator lamp flashes, press the mode selector button to switch the mode.

If the power switch is pressed after the end of the "CLEANSING" mode, the "MICRO-PAT" mode is set automatically.

3. Select the power level ("HIGH" or "LOW") with the level selector button. Select the power level ("HIGH" or "LOW") with the level selector button according to your skin condition.
 - Every time the level selector button is pressed, the power level is switched between "HIGH" and "LOW". The vibration power does not change.

Use at the "LOW" level is recommended for people who use this product for the first time and who feel something unusual.

4. Hold the main body by hand touching the side plate. (Fig. 9)
 - Make sure to hold the body with your hand touching either of the right and left side plates.

5. Apply the head directly to your skin and move it slowly. (Fig. 6)
 - The mode indicator lamp will light, and the product will start to vibrate.
 - Use this product within approximate 5 minutes.
 After 5 minutes have elapsed, the product will stop automatically.

CAUTION

The device will not vibrate unless you apply the head to your skin while touching the side plate. (Skin sensing function)
 If the head is not vibrating, the product is not working.
 If you leave the device as it is for about 30 seconds without vibration, the power will turn off automatically. (Auto-off function)

- Lightly apply and slowly move as if sliding from the center to the outside of your face.

CAUTION Do not keep the head on the same part of your skin, but move it around. (Cause of skin trouble or poor health condition)

- Do not use this device on the same part of your skin for more than 10 seconds.
- If you feel slight current stimulation (a tingling sensation) during use, press the level selector button to set the "LOW" level.

6. An automatic stop is made when around 5 minutes have elapsed. (Auto-off function)
 - You will know the automatic stop with the ending sound.
7. After use, make sure to pull out the power adaptor from the outlet. (Fig.7)
 For continued use of the next mode, keep the power adaptor connected.
8. Wipe the toner off the head and main body with tissue, etc.

Step 3 "MOISTURIZE"

This mode is to give moisturizing care by making the moisturizing ingredients of the toner permeate your skin.

1. Attach a new cotton puff to the main body, and fix it with the ring.
 - For attaching a cotton puff, see "Preparation for use."
2. Saturate the cotton puff with your regular toner.
 - Do not use a milky lotion or gel. Saturate the cotton puff with the toner until it seeps through the cotton puff to the rear side.
3. Check that the power adaptor is connected to the main body, and press the power switch. (Fig.10)
 - The "LOW" level indicator lamp will light, and the "MOISTURIZE" mode indicator lamp will flash. If another mode indicator lamp flashes, press the mode selector button to switch the mode.

If the power switch is pressed after the end of the "MICRO-PAT" mode, the "MOISTURIZE" mode is set automatically.

4. Select the power level ("HIGH" or "LOW") with the level selector button. Select the power level ("HIGH" or "LOW") with the level selector button according to your skin condition.
 - Every time the level selector button is pressed, the potential power level is switched between "HIGH" and "LOW". The vibration power does not change.

Use at the "LOW" level is recommended for people who use this product for the first time and who feel something unusual.

5. Hold the main body by hand touching the side plate. (Fig. 5)
 - Make sure to hold the body with your hand touching either of the right and left side plates.
6. Apply the cotton puff to your skin and move the product slowly. (Fig. 6)
 - The mode indicator lamp will light, and the product will start to vibrate.
 - Use this product within approximate 5 minutes.
 After 5 minutes have elapsed, the product will stop automatically.

CAUTION

The device will not vibrate unless you apply the cotton puff to your skin while touching the side plate. (Skin sensing function)
 If the head is not vibrating, the product is not working.
 If you leave the device as it is for about 30 seconds without vibration, the power will turn off automatically. (Auto-off function)

- Lightly apply and slowly move as if sliding from the center to the outside of your face.

CAUTION Do not keep the cotton puff on the same part of your skin, but move it around. (Cause of skin trouble or poor health condition)

- If the moisture in the cotton puff becomes insufficient or if it is hard to slide, add more toner.
 - Do not use this device on the same part of your skin for more than 10 seconds.
 - If you feel slight current stimulation (a tingling sensation) during use, press the level selector button to set the "LOW" level.
 If you still feel a "tingling sensation," add more toner.
7. An automatic stop is made when around 5 minutes have elapsed. (Auto-off function)
 - You will know the automatic stop with the ending sound.
 8. After use, make sure to pull out the power adaptor from the outlet. (Fig. 7)
 For continued use of the next mode, keep the power adaptor connected.
 9. Remove the ring and the cotton puff.
 - Wipe off the surface of the head with tissue, etc.
 - Do not use the cotton puff again once it is used because its reuse will cause skin trouble.

Step 4 "COOL-AESTHETIC"

This mode tightens the pores of your skin by a fresh cooling effect as the head section surface is approximately 10°C lower than the ambient temperature.

WARNING Use at an ambient temperature in the range of 10 to 40°C. (Cause of skin damage as the head is too cold if the device is used at a lower temperature)

1. Remove the cotton puff, and wipe the toner off the head section.
2. Check that the power adaptor is connected to the main body, and press the power switch. (Fig.11)
 - The level indicator lamp will go out, and the "COOL-AESTHETIC" mode indicator lamp will flash. If another mode indicator lamp flashes, press the mode selector button to switch the mode.

If the power switch is pressed after the end of the "MOISTURIZE" mode, the "COOL-AESTHETIC" mode is set automatically.

3. Wait until the head cools down.

- Preparation for cooling takes approximately 90 seconds. During preparation, the fan rotates, and the "COOL-AESTHETIC" mode indicator lamp flashes.
- When the device is ready to be used, the cooling fan stops, and a ready sound beeps, and the mode indicator lamp lights.

CAUTION

- The cooling fan may not stop even when the device is ready to be used, but this is to maintain the performance and is not a device failure.
- The head does not vibrate in the "COOL-AESTHETIC" mode.

4. While the cooling fan is operating, hold the device so as not to block the suction and exhaust holes. (Fig.12) ㉔ Suction hole ㉕ Exhaust hole

- Be careful not to block the suction and exhaust holes in the main body by hand because if the holes are blocked, the head will not cool down sufficiently.

5. Lightly apply the head directly to your skin, move and pat it slowly following the pattern of Fig.13.

- The mode indicator lamp switches to light. The head does not vibrate.
- The usable time is within approximate 90 seconds. When 90 seconds have elapsed, an automatic stop is made.
- Using at the same area for about one second, repeat "applying and releasing" the head from the center to the outside of your face.

CAUTION

- Do not slide the head as it is applied to your skin.
- Do not stop the head at the same area of your skin. The cooling effect of the head weakens.

6. When 90 seconds have elapsed, an automatic stop is made. (Auto-off function)

- The ending sound (beep) will indicate automatic stop.
- After the end of "COOL-AESTHETIC," the main body and head are slightly hot, but this is not an abnormality.

7. After use, make sure to pull out the power adaptor from the outlet. (Fig.7)

- Wipe toner, etc., if adhering to the head or main body from your skin, with tissue, etc.
- During use, the power adaptor may become slightly hot, but this is not a device failure.

CAUTION

If all the mode indicator lamps and level indicator lamps begin to flash and a warning sound beeps during use, it is a signal to indicate a device failure. Pull out the power adaptor from the outlet, and request your distributor to make repairs.

8. After use, give regular skin care with emulsion, cream, etc.

Care for the main body

- Wipe the toner or milky lotion off the main body and ring with tissue, etc. When dirt is difficult to remove, wipe it with a cloth moistened with water and firmly squeezed.

CAUTION	Do not wash or splash water on it. (Cause of failure)
	Make sure to turn OFF the switch. (Improperly turning on the switch may cause irritation or malfunction.)

- ★ **CAUTION** Do not use alcohol, benzine, or thinner. (Cause of deformation or discoloration)

How to preserve

Preserve in an area with little humidity and dust.

WARNING	Do not keep this product in places splashed with water or in damp places, such as in a bathroom. (Cause of short circuit or failure)
	Do not wind the cord around the power adaptor or main body. (Cause of damage to the cord)

- ★ **CAUTION** Do not keep this product in areas exposed to direct sunlight (inside of a vehicle in summer, etc.) or near heaters. (Cause of deformation or discoloration)

Specifications

Power supply	Power adaptor (KH-186B)	100V - 240V AC, 50-60Hz
Rated time		Around 18 minutes at the maximum
Accessories		Rings (two)

Disposal of the product

WARNING Do not disassemble except when disposing of the product.

This product includes recyclable parts.

For disposal of the parts, observe your local laws and regulations etc.

Service Center	TEL NO.
HITACHI (HONG KONG) LTD.	27535386

- Telephone numbers are subject to change without notice.

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